

ONETOUCH

Verio Reflect™

Blood Glucose Monitoring System



Owner's Booklet

Instructions for Use

ONETOUCH
Verio Reflect™
Blood Glucose Monitoring System
Owner's Booklet

Thanks for choosing OneTouch®!

The OneTouch Verio Reflect™ Blood Glucose Monitoring System is one of the latest product innovations from OneTouch®.

Your OneTouch Verio Reflect™ Meter is designed to connect (sync) with a variety of compatible wireless devices running software applications. By connecting your meter with the OneTouch Reveal® App you can review and graph your results, set Daily Test Goals, see Blood Sugar Mentor™ Messages and help identify patterns. Meter results and other tracked data are sent to the compatible devices either through BLUETOOTH® (wireless) or USB cable connection. Visit www.OnetouchMEA.com for more information about the OneTouch Reveal® App.

Every OneTouch® Meter is designed to help you test your blood glucose and manage your diabetes.

This Owner's Booklet offers a complete explanation of how to use your new meter and testing supplies. It reviews the do's and don'ts of testing your blood glucose level. Please keep your Owner's Booklet in a safe place; you may want to refer to it in the future.

We hope OneTouch® products and services will continue to be a part of your life.

Meter symbols and icons

Range Indicator icons and Colour Bar



Result below target range



Result at low end of target range



Result near middle of target range



Result in middle of target range



Result near middle of target range



Result at high end of target range



Result above target range



Result in range (text)

Meal and Event tags



Before Meal tag



After Meal tag



Event tag



Carbs Event tag



Stress Event tag



Illness Event tag



Medication Event tag



Exercise Event tag

Patterns, Test Tracker, and Awards



Low Pattern Message (blue icon)



High Pattern Message (red icon)



Daily Test Goal met



Daily Test Goal not met



Gold award



Silver award



Bronze award

Main Menu and Settings Menu icons



Summary



Results Log



Trend 90



Settings Menu



Range Indicator icons as emoji



Range Indicator icons as text

Other meter icons



BLUETOOTH® feature on



Syncing



Control solution result



Low battery



Screen navigation progress dots

Other symbols and icons

-  Cautions and Warnings: Refer to the Owner's Booklet and inserts that came with your system for safety-related information.
-  Direct current
-  Consult Instructions for Use
-  Manufacturer
-  Lot Number
-  Serial Number
-  Storage Temperature Limits
-  In Vitro Diagnostic Device
-  Do Not Re-use
-  Sterilised by irradiation
-  Not for general waste
-  Use By Date
-  Contains sufficient for n tests
-  Underwriters Laboratories certification

Before you begin

Before using this product to test your blood glucose, carefully read this Owner's Booklet, and the inserts that come with the OneTouch Verio® Test Strips, OneTouch Verio™ Control Solution and the OneTouch® Delica™ Plus Lancing Device.

⚠ CAUTION:

Healthcare professionals must also read "Multi-patient use testing". See page 142.

IMPORTANT SAFETY INSTRUCTIONS:

- After use and exposure to blood, all parts of this kit are considered biohazardous. A used kit may potentially transmit infectious diseases even after you have performed cleaning and disinfection.

Intended use

The OneTouch Verio Reflect™ Blood Glucose Monitoring System is intended for self-testing outside the body (*in vitro* diagnostic use) by people with diabetes mellitus for the quantitative measurement of glucose (sugar) in fresh capillary whole blood samples drawn from the fingertips as an aid in the management of diabetes.

The OneTouch Verio Reflect™ Meter analyses patterns, tracks trends and provides guidance (educational, behavioral and motivational messages) to help aid in the understanding and management of glucose levels and the detection of excursions above or below a desired range. The meter is also designed to wirelessly communicate with compatible diabetes applications for the capture, transfer and analysis of blood glucose data and information, to help support diabetes self-management.

The OneTouch Verio Reflect™ Blood Glucose Monitoring System should not be used as a substitute for healthcare professional advice. It also should not be used for the diagnosis or screening of diabetes or for neonatal use.

When the OneTouch Verio Reflect™ Meter is set to "Basic mode", it can be used for multi-patient testing only by healthcare professionals as an aid to monitor the effectiveness of diabetes control programmes. The system is not intended for use in critical care environments.

Test principle

Glucose in the blood sample mixes with the enzyme FAD-GDH (see page 162) in the test strip and a small electric current is produced. The strength of this current calculates your blood glucose level. It then displays the blood glucose result and stores it in the meter memory.

Use only OneTouch Verio™ Control Solution and Test Strips with the OneTouch Verio Reflect™ Meter. Use of OneTouch Verio® Test Strips with meters for which they are not intended may yield inaccurate results.

BLUETOOTH® wireless technology

BLUETOOTH® wireless technology is used by some smartphones and many other devices. Your OneTouch Verio Reflect™ Meter uses BLUETOOTH® wireless technology to pair and to send your glucose results to compatible wireless devices.

The OneTouch Verio Reflect™ Meter is designed to work with the OneTouch Reveal® Mobile App and other diabetes applications.

NOTE: Some diabetes management apps, including the OneTouch Reveal® Mobile App, may not be available in your country. Visit www.OneTouchMEA.com to learn if the OneTouch Reveal® Mobile App is available in your country.

Visit www.OneTouchMEA.com for information on which wireless devices are compatible with your OneTouch Verio Reflect™ Meter, and where/how to download the software application.

When using the OneTouch Verio Reflect™ System, we suggest you pair your OneTouch Verio Reflect™ Meter with a compatible wireless device and track your results. See page 110 for pairing instructions.

Your meter is subject to and complies with applicable Worldwide Radio regulations and guidelines. Generally, these rules state two conditions specific to the operation of the device:

1. This device may not cause harmful interference.
2. This device must accept any interference received, including interference that may cause undesirable operation.

These guidelines help ensure that your meter will not affect the operation of other nearby electronic devices. Additionally, other electronic devices should not affect the use of your meter.

If you experience meter interference problems, try moving your meter away from the source of the interference. You can also move the electronic device or its antenna to another location to solve the problem.

⚠WARNING: In locations where cell phone use is not permitted, such as hospitals, some healthcare professional offices and airplanes, you should turn the BLUETOOTH® feature off. See page 106 for more information.

Do Not leave your meter unattended or in a public place so that others may change your personal parameters or make physical connections to it.

The meter has been designed defensively to ward off malicious intent and has been tested accordingly.

BLUETOOTH® trademark

The BLUETOOTH® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by LifeScan Scotland Ltd. is under licence. Other trademarks and trade names are those of their respective owners.

Table of Contents

1	Getting to know your system.....	16
2	Setting up your system.....	26
3	Taking a test.....	40
	Testing your blood glucose.....	40
	Testing with control solution.....	66
4	Tagging results	76
5	Understanding Messages	82
6	Summary, Results Log and Trend 90	94
7	Syncing your meter	106
8	Adjusting meter settings.....	118
9	Caring for your system.....	132
10	Batteries.....	138
11	Multi-patient use testing	142
12	Troubleshooting	148
13	System information	158
14	Index	172

1 Getting to know your system

Your OneTouch Verio Reflect™ Blood Glucose Monitoring System

Included with your kit:



OneTouch Verio Reflect™
Meter (2 CR2032 lithium
coin cell batteries
included)



Lancing device



Lancets

NOTE: If any item is missing or defective in your kit, contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.

NOTE: If another type of lancing device was included, see the separate instructions for that lancing device.

1 Getting to know your system

Available separately:

Items pictured below are required, but may not be included in your kit:

They are sold separately. Refer to your meter carton for a list of included items.



OneTouch Verio™ Mid
Control Solution*



OneTouch Verio®
Test Strips*

*OneTouch Verio™ Control Solution and Test Strips are available separately. For availability of test strips and control solution, contact Customer Service or ask your pharmacist or healthcare professional.

⚠WARNING: Keep the meter and testing supplies away from young children. Small items such as the battery door, batteries, test strips, lancets, protective covers on the lancets, and control solution vial cap are choking hazards. **Do Not** ingest or swallow any items.

1 Getting to know your system

Getting to know your OneTouch Verio Reflect™ Blood Glucose Monitoring System

Meter



1	Test strip port Insert test strip here to turn meter on
2	Battery icon Appears only when batteries are low
3	mg/dL is the pre-set unit of measure and cannot be changed*
4	Micro USB data port Connect to download to a computer
5	OK button <ul style="list-style-type: none">• Turns meter on/off (press and hold)• Confirms menu selections (press and release)
6	Back button <ul style="list-style-type: none">• Go back to the previous screen (press and release)• May be used to change a Meal tag or add an Event tag
7	Up and down buttons <ul style="list-style-type: none">• Highlight selection (press and release)• Increase scrolling speed (press and hold)
8	Colour Bar
9	Range Indicator icon
10	Test Tracker icons

1 Getting to know your system

⚠WARNING: *Confirm that the unit of measure mg/dL is displayed. If your display shows mmol/L rather than mg/dL, stop using the meter and contact Customer Service.

Getting to know your OneTouch Verio® Test Strip

Test strip



11	Channel to apply sample
12	Silver prongs Insert into test strip port (1)

Features overview


The OneTouch Verio Reflect™ Meter is equipped with features designed to help you manage your blood glucose in a variety of ways. We hope you will take advantage of these features and make them part of your diabetes management routine.

Blood Sugar Mentor™ Messages

The Blood Sugar Mentor™ provides insight on how you're doing and diabetes management guidance, including educational and motivational messages based on your current and previous glucose results and other tracked data. When you receive your blood glucose result, the meter may display Mentor Tips (guidance), Pattern Messages (insight) or Awards (encouragement).

1 Getting to know your system


Examples of Blood Sugar Mentor™ Messages

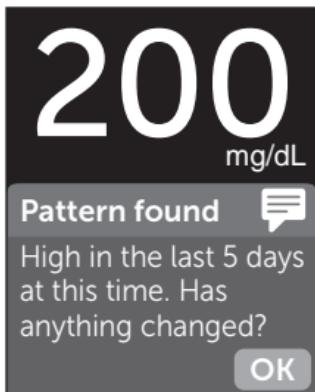

Insight



Guidance

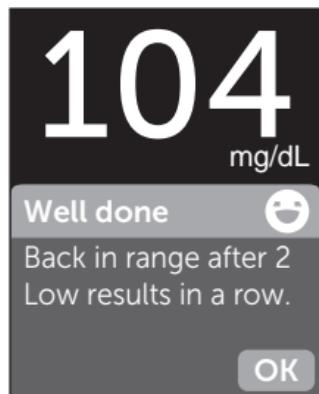

Encouragement

ColourSure™ Dynamic Range Indicator

The ColourSure™ Dynamic Range Indicator instantly lets you know if your current result is below (blue), within (green) or above (red) your range limits. It also shows when your result is near the low or high end of your target range.

Test Tracker

The Test Tracker lets you set the number of glucose tests you plan to take each day and automatically tracks progress towards achieving that goal.

Trend 90

Trend 90 lets you set your blood glucose 90 Day Average Target and tracks progress towards that goal. The graph of averages is automatically updated every two weeks and provides a view of how your overall glucose control is trending over time.

2 Setting up your system

Setting up your meter


Turn your meter on

To turn your meter on, press and hold  until the Welcome screen appears. Once the Welcome screen is displayed, release . You can also turn the meter on by inserting a test strip.

Every time you turn your meter on, a Welcome screen will appear for a few seconds. If the meter does not power on, check the batteries. See page 138.

CAUTION:

If you see any missing pixels within the Welcome screen, there may be a problem with the meter. Contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.

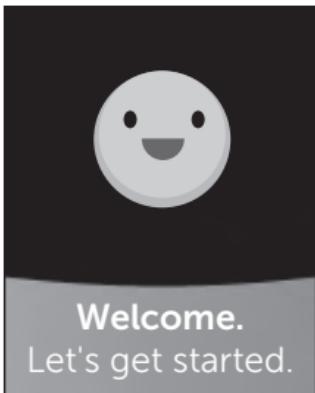

Example

NOTE: Whether you turn the meter on for the first time using  or a test strip, you will need to perform first time setup.

Use the display backlight for better visibility

The backlight comes on automatically whenever the meter is turned on. After a few seconds of no activity, the backlight will dim. Pressing any button or inserting a test strip will turn the backlight back on.

First time setup

Before using your meter for the first time, you will be prompted to set the language, time, date, range limits, Daily Test Goal, and 90 Day Average Target.

2 Setting up your system

NOTE:

- You will not be able to perform a blood glucose test until you complete first time setup.
- Be sure to talk to your healthcare professional about setting the right range limits, Daily Test Goal, and 90 Day Average Target for you. When selecting or changing these settings, you should consider factors such as your lifestyle and diabetes therapy. Never make significant changes to your diabetes care plan without consulting your healthcare professional.
- It is important to select the settings with care to ensure the correct information is stored with your test results.

The meter comes pre-set to help you with first time setup. But if you need to edit a setting, press  or  to highlight **Edit** and press . You can then press  or  to change values. After making your selection, highlight **Save** and press  to confirm your selection and proceed to the next screen.

Pressing  returns you to the previous screen.

1. Press or to highlight the language you want, then press to select

The meter will display a prompt giving the option to complete first time setup using the OneTouch Reveal® App. Visit www.OneTouchMEA.com for more information about the OneTouch Reveal® App. If you do not want this option, highlight **No** and press to continue first time setup using only the meter.



2 Setting up your system

Selecting **Yes** will prompt you to pair the meter with your compatible wireless device. See page 108. After pairing, the OneTouch Reveal® App will guide you through the process of first time setup. The remaining steps below will not be needed.

Setup

Do you want to continue setup using the OneTouch Reveal mobile app?

Yes

No




Progress dot

NOTE: The green progress dot will move from left to right as you navigate through the first time setup screens.


2. Confirm the time, then press  to save

Selecting **Edit** will allow you to change the pre-set time.

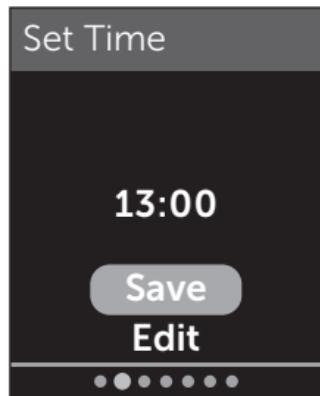

3. Confirm the date, then press  to save


Selecting **Edit** will allow you to change the pre-set date.

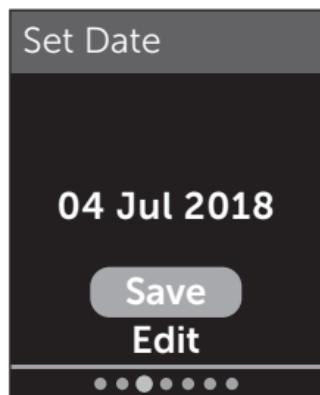

2 Setting up your system

NOTE:


- The time and date format cannot be changed.
- Once every 6 months, and every time you change the batteries, a screen will prompt you to confirm the time and date set in the meter.

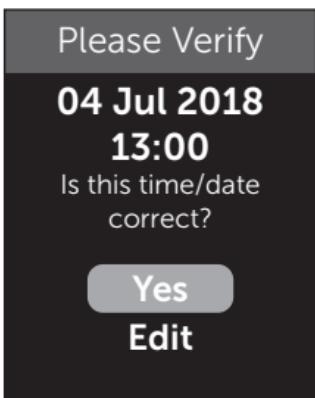

Example

Selecting **Edit** will allow you to change the time and date. Once correct, select **Done**.

4. Confirm the Before Meal range limits, then press to save

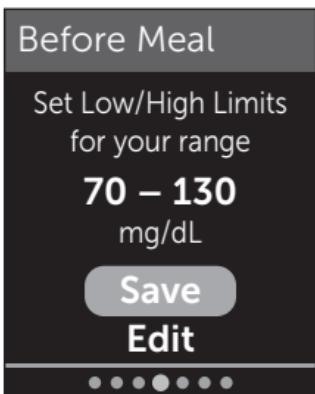

Example

Selecting **Edit** will allow you to change the pre-set Before Meal range limits.


Your meter uses the Before Meal and After Meal range limits set in your meter to tell you when a test result is within, below or above your mealtime range limits. The mealtime range limits you set during first time setup will apply to all glucose results. They are used in the ColourSure™ Dynamic Range Indicator feature, to detect Patterns and in the Summary. See page 57, page 83 and page 94 for more information.

Before Meal range limits

pre-set low limit	70 mg/dL
pre-set high limit	130 mg/dL
low limit must be set within	60 mg/dL and 110 mg/dL
high limit must be set within	90 mg/dL and 300 mg/dL

2 Setting up your system

5. Confirm the After Meal range limits, then press to save

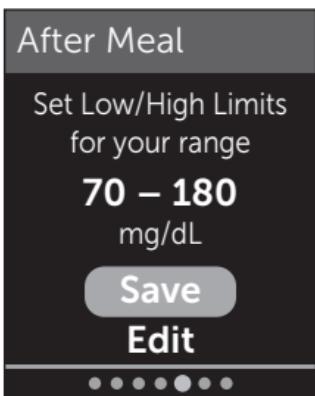

Example

Selecting **Edit** will allow you to change the pre-set After Meal range limits.

After Meal range limits


pre-set low limit	70 mg/dL
pre-set high limit	180 mg/dL
low limit must be within	60 mg/dL and 110 mg/dL
high limit must be within	90 mg/dL and 300 mg/dL

6. Confirm the Test Tracker setting, then press to save

Selecting **Edit** will allow you to change the pre-set Daily Test Goal for the Test Tracker.

The Test Tracker lets you set the number of glucose tests you plan to take each day and automatically tracks progress towards achieving that goal.

The meter is pre-set with a Daily Test Goal of four tests per day. The maximum Daily Test Goal is eight. See page 126 for more information.

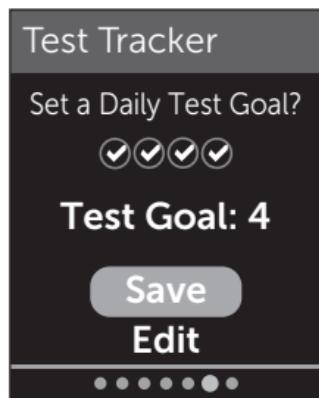


2 Setting up your system

7. Confirm the 90 Day Average Target, then press to save


Selecting **Edit** will allow you to change the pre-set 90 Day Average Target.

Trend 90 lets you set your blood glucose 90 Day Average Target and tracks progress towards that goal. The graph of averages is automatically updated every two weeks and provides a view of how your overall glucose control is trending over time.

The meter is pre-set with a 90 Day Average Target of 155 mg/dL. The minimum 90 Day Average Target is 99 mg/dL and the maximum is 265 mg/dL. See page 103 for more information.

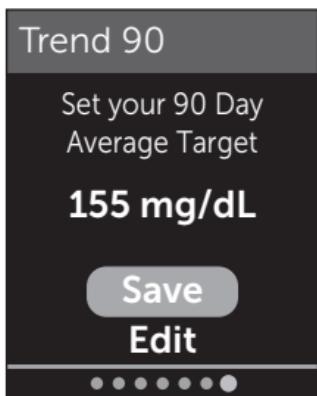

Setup complete appears on the screen. Your meter is ready for use.

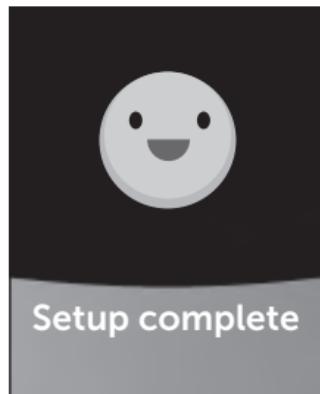

If the meter was turned on by inserting a test strip, the **Apply Blood** screen will appear. To perform a blood glucose test, see page 40.

NOTE: After first time setup, you can adjust your settings at any time. See page 118 for more information.

Understanding Basic Mode

The meter comes with many features to help track diabetes management progress. Basic Mode will turn off these features if you prefer not to use the meter in this way. Healthcare professionals should select Basic Mode to disable features that are inappropriate for a clinical multi-patient setting.

2 Setting up your system

When Basic Mode is selected, the Test Tracker, Meal and Event tagging, Before and After Meal averages, Trend 90, Pattern Messages, Mentor Tips and Awards will be turned off. In Basic Mode, range limits are general and are not defined by mealtime. Range Indicator icons appear as text (for example, ) rather than as emoji (for example, ).

For more information about turning Basic Mode on or off see page 125.

Turning the meter off

There are three ways to turn your meter off:

- Press and hold  for several seconds until the meter turns off.
- Remove the test strip.
- Your meter will turn off by itself if left alone for two minutes.

This page left blank intentionally.

3 Taking a test

Testing your blood glucose

Preparing for a blood glucose test

NOTE: Many people find it helpful to practise testing with control solution before testing with blood for the first time. See page 68.

Have these things ready when you test:

OneTouch Verio Reflect™ Meter

OneTouch Verio® Test Strips

Lancing device

Sterile lancets

NOTE:

- Use only OneTouch Verio® Test Strips.
- Unlike some blood glucose meters, no separate step to code your OneTouch Verio Reflect™ System is required.
- Testing must be done within operating temperature range (6-44°C).
- Make sure your meter and test strips are about the same temperature before you test.

- Keep test strips in a cool, dry place between 5°C and 30°C.
- **Do Not** test if there is condensation (water build-up) on your meter. Move your meter and test strips to a cool, dry spot and wait for the meter surface to dry before testing.
- **Do Not** open the test strip vial until you are ready to remove a test strip and perform a test. Use the test strip **immediately** after removing it from the vial.
- Tightly close the cap on the vial immediately after use to avoid contamination and damage.
- Store unused test strips only in their original vial.
- **Do Not** return the used test strip to the vial after performing a test.
- **Do Not** re-use a test strip that had blood, control solution, or any contaminants applied to it. Test strips are for single use only.
- With clean, dry hands, you may touch the test strip anywhere on its surface. **Do Not** bend, cut or modify the test strip in any way.

3 Taking a test

- When you first open a vial of test strips, record the discard date on the label. Refer to the test strip insert or vial label for instructions on determining the discard date.
- Comparing your blood glucose test results taken with this meter to your results taken from a different meter is not recommended. Results may differ between meters and are not a useful measure of whether your meter is working properly. To check your meter accuracy, you should periodically compare your meter results to those obtained from a lab. See page 158 for more information.

IMPORTANT: If another person assists you with testing, the meter should always be cleaned prior to use by that person. See page 132.


⚠ CAUTION:

- The OneTouch Verio Reflect™ Blood Glucose Monitoring System should not be used for patients within 24 hours of receiving a D-xylose absorption test as it may cause inaccurately high results.

- **Do Not** use the OneTouch Verio™ Family of Meters when PAM (Pralidoxime) is known or suspected to be in the patient's whole blood sample.
- **Do Not** use your test strips if your vial is damaged or left open. This could lead to error messages or inaccurate results. Contact Customer Service immediately if the test strip vial is damaged. Contact information for Customer Service is listed at the end of this Owner's Booklet.
- If you cannot test due to a problem with your testing supplies, contact your healthcare professional. Failure to test could delay treatment decisions and lead to a serious medical condition.
- The test strip vial contains drying agents that are harmful if inhaled or swallowed and may cause skin or eye irritation.
- **Do Not** use test strips after the expiry date (printed on the vial) or the discard date, whichever comes first, or your results may be inaccurate.

3 Taking a test

Getting to know your OneTouch® Delica™ Plus Lancing Device

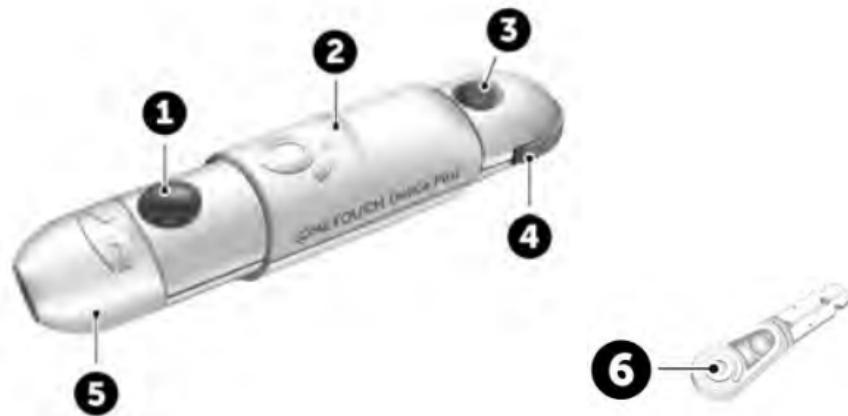

1	Release button
2	Slider control
3	Depth indicator
4	Depth wheel
5	Lancing device cap
6	Protective cover

NOTE:

- The OneTouch® Delica™ Plus Lancing Device uses OneTouch® Delica™ or OneTouch® Delica™ Plus Lancets.
- If another type of lancing device was included, see the separate instructions for that lancing device.
- The OneTouch Verio Reflect™ Blood Glucose Monitoring System has not been evaluated for Alternate Site Testing (AST). Use only fingertips when testing with the system.
- The OneTouch® Delica™ Plus Lancing System does not include the materials needed to perform Alternate Site Testing (AST). The OneTouch® Delica™ Plus Lancing System should not be used on the forearm or palm with the OneTouch Verio Reflect™ Blood Glucose Monitoring System.

Lancing precautions

⚠ CAUTION:


To reduce the chance of infection and disease spread by blood:

- Make sure to wash the sample site with soap and warm water, rinse and dry before sampling.
- The lancing device is intended for a single user. Never share a lancet or lancing device with anyone.
- Always use a new, sterile lancet each time you test.
- Always keep your meter and lancing device clean. (See page 132.)
- After use and exposure to blood, all parts of this kit are considered biohazardous. A used kit may transmit infectious diseases even after you have performed cleaning and disinfection.
- **Do Not** use lancets after the expiry date printed on the lancet packaging.

Preparing your lancing device


1. Remove the lancing device cap

Remove the cap by rotating it and then pulling it straight off the device.

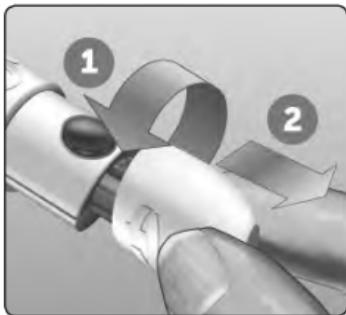

2. Insert a sterile lancet into the lancing device

Align the lancet as shown here, so that the lancet fits into the lancet holder. Push the lancet into the device until it snaps into place and is fully seated in the holder.

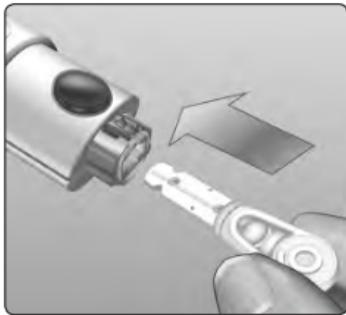

3 Taking a test

Twist the protective cover one full turn until it separates from the lancet. **Save the protective cover for lancet removal and disposal.**



See page 63.



3. Replace the lancing device cap


Place the cap back onto the device; turn the cap or push the cap straight on to secure it.

Ensure the cap is aligned as shown in the image.

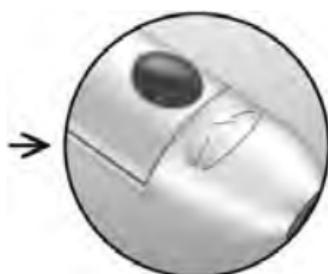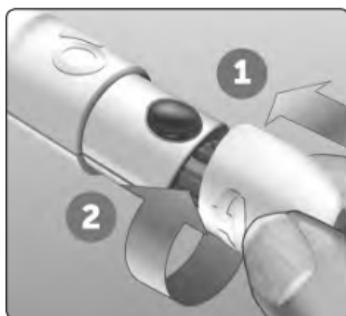

4. Adjust the depth setting

The device has 13 puncture depth settings (each dot shown between numbers 1 to 7 on the Depth Wheel indicates an additional available depth setting). Adjust the depth by turning the depth wheel. Smaller numbers are for a shallower puncture and larger numbers for a deeper puncture.

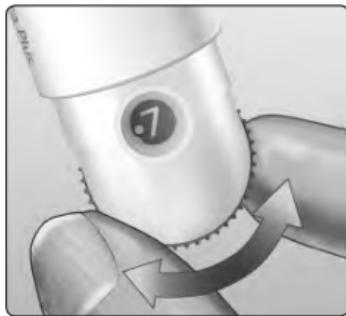

NOTE: Try a shallower setting first and increase the depth until you find the one deep enough to get a blood sample of the proper size.

5. Cock the lancing device

Pull the slider back until it clicks. If it does not click, it may already have been cocked when you inserted the lancet.



3 Taking a test

Preparing the meter

1. Insert a test strip to turn the meter on

Insert a test strip into the test strip port with the gold side of the test strip and the two silver prongs facing you.


Silver prongs



Test strip port

NOTE: No separate step to code the meter is required.

The **Apply Blood** screen appears on the display. You can now apply your blood sample to the test strip.

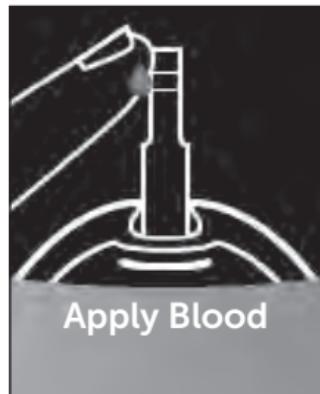

Sampling blood from the fingertip


Choose a different puncture site each time you test. Repeated punctures in the same spot may cause soreness and calluses.

Before testing, wash your hands and the sample site with warm, soapy water. Rinse and dry completely. Contaminants on the skin may affect results.

3 Taking a test


1. Puncture your finger

Hold the lancing device firmly against the side of your finger. Press the release button. Remove the lancing device from your finger.

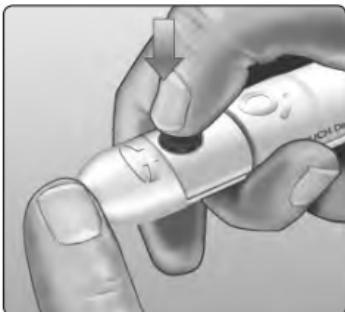

2. Get a round drop of blood

Gently squeeze and/or massage your fingertip until a round drop of blood forms.

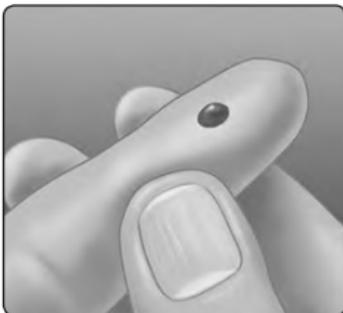

NOTE: If the blood smears or runs, **Do Not** use that sample. Dry the area and gently squeeze another drop of blood or puncture a new site.



Applying blood and reading results

Applying the sample


NOTE: Healthcare professionals must also read additional instructions in Multi-patient use testing. See page 142.

1. Apply the sample to the test strip

You can apply blood to either side of the test strip.


Apply your sample to the opening of the channel.

Be sure to apply your sample immediately after you get a drop of blood.

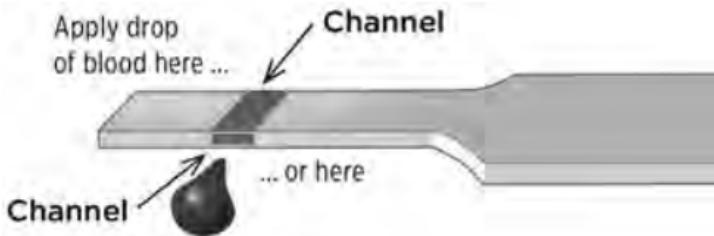


3 Taking a test

Holding the meter at a slight angle, guide the channel to the blood drop.

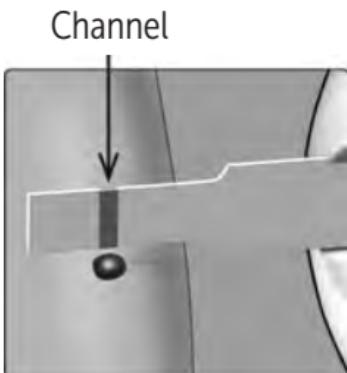

When it touches your sample, the test strip wicks blood into the channel.

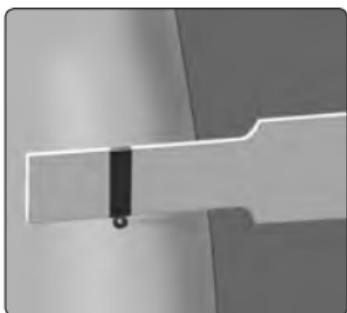


2. Wait for the channel to fill completely

The blood drop will be drawn into the narrow channel. The channel should fill completely.

The channel turns red and a result will be displayed after approximately 5 seconds.

Blood should **not** be applied on the top of the test strip or to the top edge of the test strip.

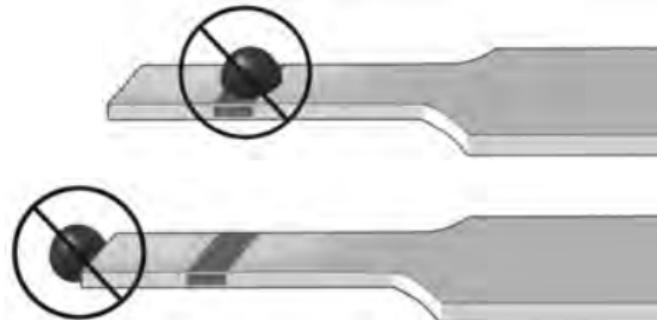

3 Taking a test

NOTE:

- **Do Not** smear or scrape the sample with the test strip.
- **Do Not** press the test strip too firmly against the puncture site or the channel may be blocked from filling properly.
- **Do Not** apply more blood to the test strip after you have moved the drop of blood away.
- **Do Not** move the test strip in the meter during a test or you may get an error message or the meter may turn off.
- **Do Not** remove the test strip until the result is displayed or the meter will turn off.

Viewing your result

ColourSure™ Dynamic Range Indicator

The OneTouch Verio Reflect™ Meter comes with a ColourSure™ Dynamic Range Indicator feature that instantly lets you know if your current result is below (blue), within (green) or above (red) your range limits. It also shows when your result is near the low or high end of your target range. It does this by displaying your result with a Range Indicator icon pointing to a segment on the Colour Bar based on the ranges you have set in the meter. The Range Indicator icon may be either an emoji (for example, 😊) or text (for example **In Range**). Use the Range Indicator icon and the segmented Colour Bar together to see where your result falls within your range limits.

3 Taking a test



Low
(Blue)

Example - below
range result



In Range
(Green)

Example - low end of
range result



High
(Red)

Example - above
range result

If the Test Tracker is turned on, progress toward the Daily Test Goal will be indicated above the result. If the Test Tracker is turned off, the date and time will be shown above the result. The date and time are always recorded with each result and can be viewed in the Results Log. See page 100.


⚠WARNING: Confirm that the unit of measure mg/dL is displayed. If your display shows mmol/L rather than mg/dL, stop using the meter and contact Customer Service.

⚠CAUTION:

Do Not make immediate treatment decisions based on the ColourSure™ Dynamic Range Indicator feature. Treatment decisions should be based on the numerical result and healthcare professional's recommendation and not solely on where your result falls within your range limits.

3 Taking a test

After the result screen appears, a Pattern Message, Mentor Tip or Award may be displayed. The Message can be dismissed by pressing . If you would like to bring the Message back, press  again. See page 82 for more information about Messages.

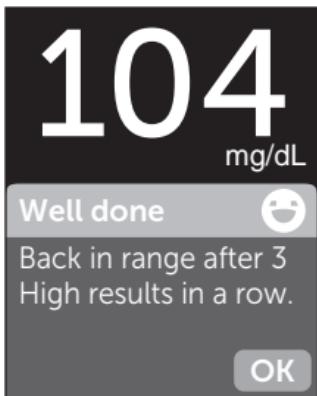

Example


A Before Meal tag will automatically be applied to your result. When you change a Before Meal tag to an After Meal tag, the Range Indicator icon may point to a different segment on the Colour Bar. A new Message may appear. See page 76 for information on tagging results.

Interpreting unexpected results

Refer to the following cautions when your results are higher or lower than what you expect.

⚠️CAUTION:**Low blood glucose results**

If your blood glucose result is below 70 mg/dL or is shown as **EXTREME LOW GLUCOSE**, (meaning the result is less than 20 mg/dL), it may mean hypoglycaemia (low blood glucose). This may require immediate treatment according to your healthcare professional's recommendations. Although this result could be due to a test error, it is safer to treat first, then do another test.

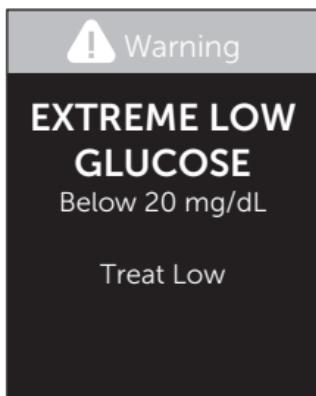**⚠️CAUTION:****Dehydration and low blood glucose results**


You may get false low blood glucose results if you are severely dehydrated. If you think you are severely dehydrated, contact your healthcare professional immediately.

3 Taking a test

⚠️ CAUTION:

High blood glucose results

If your blood glucose result is above 180 mg/dL, it may mean hyperglycaemia (high blood glucose) and you should consider re-testing. Talk to your healthcare professional if you are concerned about hyperglycaemia.

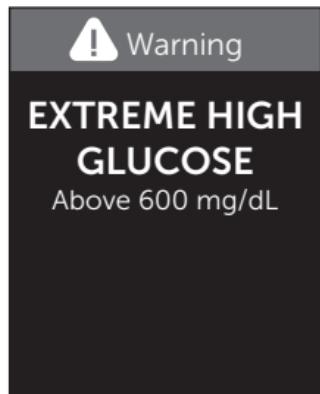

EXTREME HIGH GLUCOSE is displayed when your blood glucose result is over 600 mg/dL. You may have severe hyperglycaemia (very high blood glucose). Re-test your blood glucose level. If the result is **EXTREME HIGH GLUCOSE** again, this indicates a severe problem with your blood glucose control. Obtain and follow instructions from your healthcare professional immediately.

⚠CAUTION:**Repeated unexpected blood glucose results**

- If you continue to get unexpected results, check your system with control solution. See page 68.
- If you are experiencing symptoms that are not consistent with your blood glucose results and you have followed all instructions in this Owner's Booklet, call your healthcare professional. Never ignore symptoms or make significant changes to your diabetes management programme without speaking to your healthcare professional.

Unusual red blood cell count

A haematocrit (percentage of your blood that is red blood cells) that is either very high (above 60%) or very low (below 20%) can cause false results.


Removing the used lancet

NOTE: This lancing device has an ejection feature, so you do not have to pull out the used lancet.

3 Taking a test


1. Remove the lancing device cap

Remove the cap by rotating it and then pulling it straight off the device.

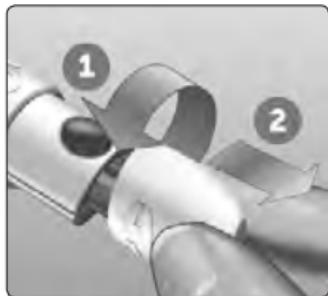

2. Cover the exposed lancet tip

Before removing the lancet, place the lancet protective cover on a hard surface then push the lancet tip into the flat side of the disk.

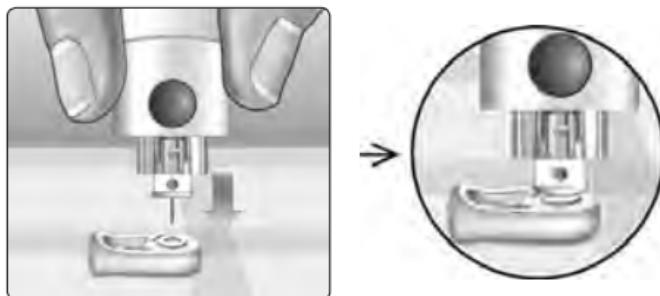

3. Eject the lancet



Holding the lancing device directed downwards, push the slider forward until the lancet comes out of the lancing device. If the lancet fails to eject properly, cock the device then push the slider control forward until the lancet comes out.



4. Replace the lancing device cap

Place the cap back onto the device; turn the cap or push the cap straight in to secure it.

Ensure the cap is aligned as shown in the image.

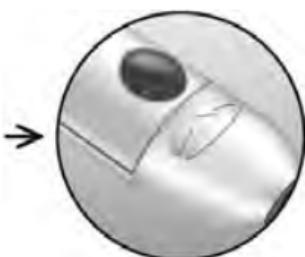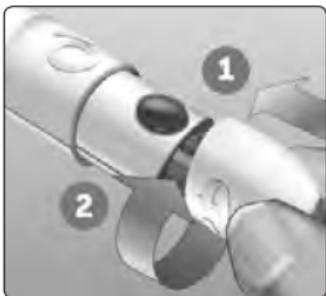

3 Taking a test

It is important to use a new lancet each time you obtain a blood sample. **Do Not** leave a lancet in the lancing device. This will help prevent infection and sore fingertips.

Disposing of the used lancet and test strip

Discard the used lancet carefully after each use to avoid unintended lancet stick injuries. Used lancets and test strips may be considered biohazardous waste in your area. Be sure to follow your healthcare professional's recommendations or local regulations for proper disposal.

Wash hands thoroughly with soap and water after handling the meter, test strips, lancing device and cap.

Testing with control solution

Control solution testing precautions

OneTouch Verio™ Control Solution is used to check that the meter and test strips are working together properly and that the test is performing correctly. (Control solution is available separately.)

NOTE:

- Use only OneTouch Verio™ Control Solution with your OneTouch Verio Reflect™ Meter.
- When you first open a new vial of control solution, record the discard date on the vial label. Refer to the control solution insert or vial label for instructions on determining the discard date.
- Tightly close the cap on the control solution vial immediately after use to avoid contamination or damage.

⚠ CAUTION:

- **Do Not** swallow or ingest control solution.
- **Do Not** apply control solution to the skin, eyes, ears or nose as it may cause irritation.
- **Do Not** use control solution after the expiry date (printed on the vial label) or the discard date, whichever comes first, or your results may be inaccurate.

3 Taking a test

When to do a control solution test

- When you open a new vial of test strips.
- If you suspect that the meter or test strips are not working properly.
- If you have had repeated unexpected blood glucose readings.
- If you drop or damage the meter.

Preparing your meter for a control solution test

1. Press and hold to turn the meter on


Wait for the Main Menu to appear.

2. Press or to highlight Settings, then press



3. Highlight Control Test, then press 

Wait for the **Insert Strip** screen to appear on the display.

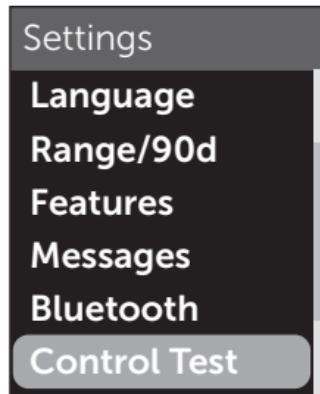

4. Insert a test strip into the test strip port



5. Wait for the Apply Control Solution screen to appear on the display



3 Taking a test

Preparing the control solution

1. Before removing the cap, shake the vial gently
2. Remove the vial cap and place it on a flat surface with the top of the cap pointing up
3. Squeeze the vial to discard the first drop



5. Squeeze a drop into the small well on the top of the cap or onto another clean, non-absorbent surface



Applying the control solution

1. Hold the meter so that the side edge of the test strip is at a slight angle to the drop of control solution



2. Touch the channel on the side of the test strip to the control solution

3. Wait for the channel to fill completely



3 Taking a test

Viewing your control solution result

After the control solution is applied, your meter will display a progress screen for approximately 5 seconds. Your result is displayed along with the date, time, unit of measure, and **C** (for control solution).



Example

Control solution results are stored in the meter and can be seen when reviewing past results on the meter.

⚠ CAUTION:

Make sure you select **Control Test** from the Settings Menu before you begin a control solution test. Follow the steps beginning on page 68. If **C** does not appear on the screen, this result will be included in your averages and your averages will change too. Repeat the test with a new test strip. If the problem persists, Contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.

Checking if the result is in range

Compare the result displayed on the meter to the range printed **on your OneTouch Verio™ Control Solution vial**. If your control solution result falls outside the expected range, repeat the test with a new test strip.



Example range
OneTouch Verio™ Mid
Control Solution Control
Range 102-138 mg/dL

⚠ CAUTION:

The range printed on the control solution vial is for control solution tests only **and is not a recommended range for your blood glucose level**.

3 Taking a test

Causes of out-of-range results

Out-of-range results may be due to:

- Not following the instructions for performing a control solution test.
- Control solution is contaminated, expired, or past its discard date.
- Test strip or test strip vial is damaged, expired, or past its discard date.
- Meter, test strips and/or control solution were not all at the same temperature when the control solution test was performed.
- A problem with the meter.
- Dirt or contamination in the small well on the top of the control solution cap.

⚠CAUTION:

If you continue to get control solution results that fall outside the range printed on the control solution vial, **Do Not** use the meter, test strips, or control solution. Contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.

Cleaning the control solution cap

Clean the top of the control solution cap with a clean, damp tissue or cloth.

4 Tagging results

Adding Meal and Event tags

Your OneTouch Verio Reflect™ Meter allows you to attach a Before Meal (🍎) or After Meal (🍏) tag to your current blood glucose result. It is important to understand the link between blood glucose levels and food. When you tag your test results with Meal tags, you add information about food to your results. Your meter uses this information in a variety of ways and displays Pattern Messages, Mentor Tips and Awards. Talk to your healthcare professional to see how Meal tags may help you manage your diabetes. Consider the following guidelines when Meal tagging.

- A Before Meal blood glucose test is taken just before the start of your meal.
- An After Meal blood glucose test is typically taken one to two hours after the start of your meal.
- A blood glucose test taken while fasting should be tagged as Before Meal.

In addition to Meal tagging, your OneTouch Verio Reflect™ Meter includes an Event tagging feature. Event tagging adds information about your health and lifestyle events to your blood glucose results. Your meter will review your past results and display Messages associated with tagged results when certain conditions are met. See page 4 for Event tag icons and their meanings.

1. View the current blood glucose result on the display

The result will automatically receive a Before Meal tag and no Event tag.

2. If a Message appears, press  to clear the Message

4 Tagging results

3. To change to an After Meal tag, press (OK) then highlight After Meal and press (OK)

NOTE: The Range Indicator icon displayed with the result uses Before and After Meal ranges. When you change a Before Meal tag to an After Meal tag, the Range Indicator icon may point to a different segment on the Colour Bar. A new Message may appear.



4. To add an Event tag to your result, press (OK)

5. Press  or  to highlight Event Tag, then press 



6. Press  or  to highlight an Event tag, then press 

You may add more than one Event tag to a result.

A  notes that an Event tag has been added.



4 Tagging results

7. When you have finished selecting Event tags, highlight Done and press 



Example

The result will be displayed with the tags. Pressing and holding  will allow you to edit the tags.

NOTE: Use care when adding tags to blood glucose results. When you leave the result screen, the tags and result are stored in the Results Log and **cannot** be changed. Incorrect Meal tagging can cause the meter to identify inaccurate Before and After Meal averages and Patterns. Incorrect Event tagging can cause the meter to display incorrect Messages.

This page left blank intentionally.

5 Understanding Messages

Blood Sugar Mentor™ Messages overview

The OneTouch Verio Reflect™ Meter reviews your past test results and may display Mentor Tips (guidance), Pattern Messages (insight) or Awards (encouragement). The meter will show you the most appropriate Message if more than one is applicable to your result.

1. **Mentor Tips** appear when your results are consistently in range or are currently trending low or high.
2. **Pattern Messages** appear when the meter identifies a pattern of glucose results that fall outside the high and low range limits you set in the meter.
3. **Awards** are earned when certain Award criteria are met, such as meeting the Daily Test Goal or tagging results.

Pressing  will dismiss Messages.

Understanding Pattern Messages

Low and High Pattern Messages appear when the meter identifies a Pattern of glucose results that fall outside the low and high mealtime range limits set in your meter.

Every time you test your blood glucose, your OneTouch Verio Reflect™ Meter searches for any new Patterns that have developed over the past 5 days. Patterns are identified by finding results below or above the low and high mealtime range limits set in your meter, and cross referencing with the time of day the tests were taken. For a set of results to form a Pattern, the results must be within the same 3-hour time period over the past 5 days.

When a Pattern is identified, a Pattern message icon (💬) appears below your result. A Pattern Message will be displayed after the result screen. Inaccurate results may cause Pattern Messages to appear.

Once a result is used in a Pattern, it will not be used again in future Pattern Messages.

5 Understanding Messages


Reviewing results that create Patterns

There are two ways you can review the individual results that combine to create a Low or High Pattern.

To view a Low or High Pattern message after a test, press .



Example - Results screen with Pattern icon

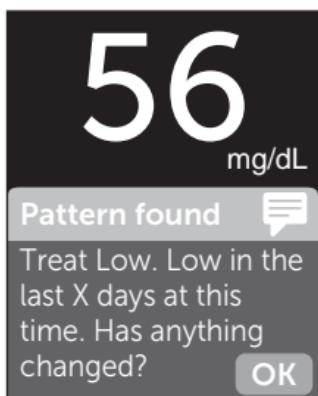

Example - Pattern Message

To view the details about the individual blood glucose results that resulted in the Pattern, press  again.

Low Pattern		
mg/dL		
11 Jun		56
00:05		
09 Jun		65
00:01		

Example - Results that created Low Pattern

Or you can access the Patterns using the Results Log screen. See page 100.

5 Understanding Messages

Low (or Hypo*) Pattern Messages appear when any 2 results over the past 5 days are:

- within the same 3-hour period
- below your low mealtime range limits

NOTE: Follow your healthcare professional's advice for treating a low blood glucose result.

*Lee-Davey, J., Alexander, S., & Raja, P. (2011, February 16). *Clinical Patterns of Low Blood Glucose Identified by a Pattern Algorithm May Predict Increased Risk of Severe Hypoglycaemia in the Following 24-Hour Period* [PDF]. Inverness, Scotland: LifeScan.

High Pattern Messages appear when any 3 results over the past 5 days are:

- within the same 3-hour period
- above your high mealtime range limits
- share the same Meal tag

Event tag Pattern Messages appear when any 3 results over the past 30 days are:

- within the same 3-hour period
- above or below your mealtime range limits
- share the same Event tag

5 Understanding Messages

NOTE: To be sure that Low and High Pattern Messages appear:

- Basic Mode must be turned off and Pattern Messages must be turned on. See page 125 and page 128.
- Update the date and time if you change time zones, or if the time changes due to Daylight saving time.
- Test your blood glucose using only one OneTouch Verio Reflect™ Meter. Using multiple meters may cause you to miss Patterns.
- Tag your blood glucose results carefully.
- Test when you feel your blood glucose may be high or low.

⚠CAUTION:

- Always use your current result for immediate treatment decisions.
- **Do Not** use Pattern Messages to make immediate or significant changes to your diabetes care plan without first consulting your healthcare professional. Always consult your healthcare professional before making significant changes to your diabetes care plan.
- **Do Not** wait for Pattern Messages to treat low or high results.

5 Understanding Messages

Understanding Mentor Tips

Mentor Tips let you know when your results are consistently in range, and display diabetes management information when results are trending low and high.

NOTE: Mentor Tips appear with current results and cannot be viewed later.

- **Unusually Low** - Current result is below the low range limit, and is at least 5 mg/dL below your lowest result over the past 14 days.
- **Morning Low** - Current result occurred in the morning and is below the low range limit.¹
- **Treat Low** - Current result is below the low range limit.^{1, 2}
- **Unusually High** - Current result is above the high range limit, and is more than 10% higher than the highest blood glucose result in the last 14 days.

- **Morning High** - Current result occurred in the morning and is above the high range limit.¹
- **Bedtime High** - Current result occurred at night and is above the high range limit.¹
- **Congratulations** - Your 90-day average is improving.
- **Trend 90 Rise** - Your 90-day average has increased.
- **Trend 90 Update** - New 90-day average is available.
- **Near Low** - Result is in range but close to low range limit.
- **Near High** - Result is in range but close to high range limit.
- **Well Done** - Result is back in range after a series of either 2 low or 3 high results.
- **Keep It Up** - More than 70% of results over the past 7 days were in range.
- **Weekly Average** - New 7-day average is shown.

5 Understanding Messages

- **Device Not Paired** - Reminds you to pair the meter with a compatible wireless device.
- **Connect Device** - Reminds you to sync the meter with the App.
- **Exercise Low** - Current result is below the low range limit and received an Exercise Event tag.^{3, 4}
- **Take a Break** - Current result is above the high range limit and received a Stress Event tag.⁴

¹Riddle, M. C., MD (Ed.). (2018). 6. Glycemic Targets: Standards of Medical Care in Diabetes 2018 [Abstract]. *Diabetes Care*, 41 (Suppl. 1), S60-61.

²American Association of Diabetes Educators. (2018). *HYPOGLYCEMIA What are the symptoms? And how should I treat it?* [Brochure]. Chicago, IL: Author.

³American Association of Diabetes Educators. (2017). *AADE7 Self-Care Behaviors, Problem Solving* [Brochure]. Chicago, IL: Author.

⁴Riddle, M. C., MD (Ed.). (2018). 6. Glycemic Targets, S44-45.

Understanding Awards

The OneTouch Verio Reflect™ Meter tracks the goals outlined below. Each time you earn bronze (bronze medal icon), silver (silver medal icon) or gold (gold medal icon) medal status, an Award will appear along with your in-range result. You can review your Awards later in the Summary. (See page 94.)

- **Test Goal Award** - Meet your Daily Test Goal a set number of times in total or in a row.
- **In-Range Award** - Achieve an in-range result a set number of times in a row.
- **Event Tag Award** - Event tag your result a set number of times in total.
- **Syncing Award** - Sync your meter to the OneTouch Reveal® App a set number of times in total.

6 Summary, Results Log and Trend 90


Viewing the Summary

Depending on settings selected the Summary may display results by range, time of day, averages, mealtime averages, Test Tracker progress and Awards. You and your healthcare professional can reference the Summary to identify overall trends in your diabetes management.

NOTE:

- **Do Not** use result averages to make immediate treatment decisions. Result averages provide information from past results. Always consult your healthcare professional before making significant changes to your diabetes care plan.
- **Do Not** allow other people to use your meter as it may affect your averages.
- The meter calculates averages based on the 7-, 14- and 30-day periods ending on the current date setting. If you change the date setting, the averages may change too.

- The Summary screens displayed will reflect the features currently selected. See page 118 for information on turning meter features on or off.
- If Basic Mode is selected, your Summary will be limited to the 30-day Summary and Average views. See page 37 for more information about Basic Mode.
- In Averages, Pattern Messages and Awards, an **EXTREME HIGH GLUCOSE** result is always counted as 600 mg/dL, and an **EXTREME LOW GLUCOSE** result is always counted as 20 mg/dL.
- The green progress dot will move from left to right as you navigate through the summary screens.

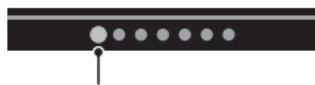


- Pressing will display the previous screen.

6 Summary, Results Log and Trend 90

1. From the Main Menu, press or to highlight Summary and press




The number of results that are low (blue), in-range (green) and high (red) are displayed and appear as bar graphs.

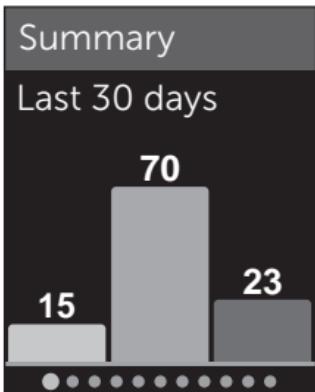

Example

2. Press  to display the Summary by time of day

The number of test results that are low, in-range and high, by time of day are displayed. The four time slots are Morning, Afternoon, Evening and Bedtime. The time frames for the four time slots cannot be changed.

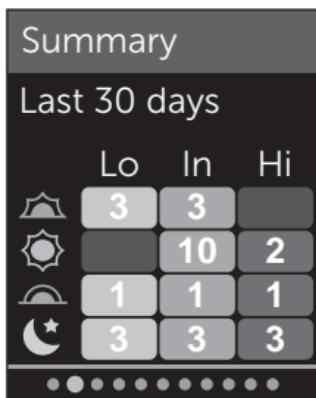

Example

Morning 	06:00 to 11:59
Afternoon 	12:00 to 17:59
Evening 	18:00 to 23:59
Bedtime 	00:00 to 05:59

NOTE: If there were no results in the past 30 days in a target range and time period, that section of the chart is left blank.

6 Summary, Results Log and Trend 90

3. Press  to display your overall averages for the past 7, 14 and 30 days

Averages	
	mg/dL
7 Days	160
14 Days	180
30 Days	200

Example

4. If Basic Mode is off, press  to display your Before Meal averages for the past 7, 14 and 30 days

Before Meal Avg.
mg/dL
7 Days  200
14 Days  ---
30 Days  ---

Example

5. If Basic Mode is off, press  to display your After Meal averages for the past 7, 14 and 30 days

After Meal Avg.	
mg/dL	
7 Days	 160
14 Days	 180
30 Days	 200
• • • • • • •	

Example

NOTE: If you do not have results in the past 7-, 14- and 30-day periods, dashes will appear in the mg/dL column.


6. If Test Tracker is on, press  to display the Daily Test Goal screens

Summary
Test Tracker
 
4 out of 8 test(s) completed today.

Example

6 Summary, Results Log and Trend 90

7. If Awards is on, press  to display the Award screens

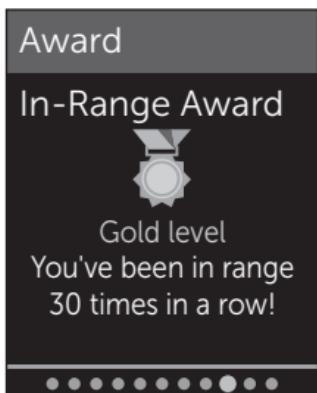

Example

Viewing the Results Log

Using the Results Log, you can review the most recent 750 results.

See page 4 for definitions of meter icons that may appear with your result.

Pressing and releasing  will allow you to view the previous screen.

1. From the Main Menu, press  or  to highlight Results Log and press 



Up to four results are displayed on the screen, starting with the most recent.

Results Log	
mg/dL	
Mon, 11 Jun	
18:05	 304
00:05	 56
Sun, 10 Jun	
19:05	 150

Example

6 Summary, Results Log and Trend 90

2. Press or to move through the list of results

NOTE: If a result is part of a Pattern, the Pattern icon will be displayed for that result instead of any tags that result may have.

Results Log	
mg/dL	
Mon, 11 Jun	
18:05	304
00:05	56
Sun, 10 Jun	
19:05	150

Example

3. Highlight a result with an icon and press

The result screen will be displayed along with the date and time the result was taken and any tag, Pattern or Range Indicator icons. If you are uncertain whether the result you are viewing is your most current result, view the date and time on the screen.



Example

If the result is part of a Pattern, indicated by  , pressing  again will display the results that led to the Pattern. See page 83 for more information about Patterns.

Low Pattern		mg/dL
11 Jun	00:05	 56
09 Jun	00:01	 65

Example

Viewing the Trend 90


The OneTouch Verio Reflect™ Meter averages results over the past 90 days to provide a view of the overall trend of blood glucose results. There must be at least 90 days of results stored in the meter for the first 90-day average to appear. A point will be graphed on the screen to represent the first 90-day average. After the first 90-day average a new point is displayed every 14 days.

6 Summary, Results Log and Trend 90


NOTE:

- Basic Mode must be turned off to use the Trend 90 feature. (See page 125.)
- The Trend 90 is an average of blood glucose results. Always use your current result for immediate treatment decisions.

**1. From the Main Menu,
press  or  to highlight
Trend 90 and press **

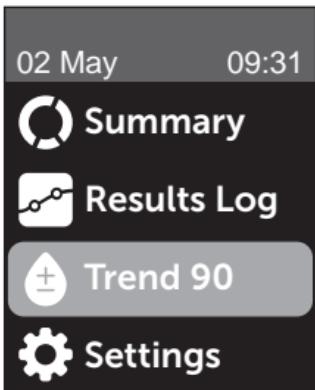

The current 90-day average will be displayed as a blood glucose value with the unit of measure. Each point on the graph represents a 90-day average automatically generated every 14 days. You can track your Trend 90 visually by following the path of points along the graph.

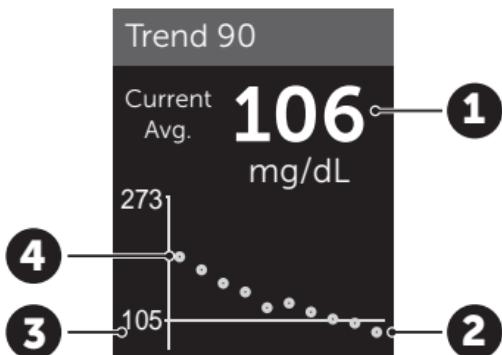

Example

1	Current 90-day average
2	Newest average point
3	90 Day Average Target set in the meter
4	Oldest average point

2. Press to return to the Main Menu

7 Syncing your meter

Connecting to a compatible wireless device


Turning the BLUETOOTH® feature on or off

In order to connect your meter with your compatible wireless device, the BLUETOOTH® feature will need to be turned on. The  symbol will appear on the meter screen when the BLUETOOTH® feature is on. When the  symbol is not present on the screen the BLUETOOTH® feature is off.

1. Start by turning your meter on using the button

Wait for the Main Menu to appear, then press  or  to highlight **Settings** and press 

2. On the Settings Menu, highlight Bluetooth and press

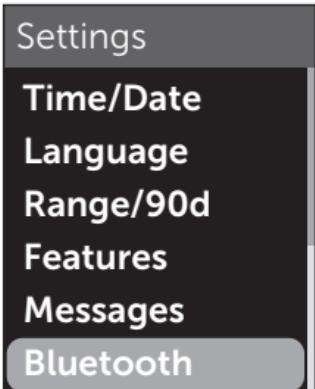

A ✓ notes if BLUETOOTH® is currently set to On or Off.



3. Press or to highlight On or Off and press

NOTE: The BLUETOOTH® feature will turn OFF during a blood glucose test.

7 Syncing your meter

Pairing overview

Pairing allows your OneTouch Verio Reflect™ Meter to communicate with compatible wireless devices. The devices must be within 8 metres of each other to pair and sync. Download the OneTouch Reveal® Mobile App from the App Store or on Google Play before pairing your meter and compatible wireless device.



Download on the
App Store



GET IT ON
Google Play

Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

Google Play and the Google Play logo are trademarks of Google LLC.

NOTE: Some diabetes management apps, including the OneTouch Reveal® Mobile App, may not be available in your country. Visit www.OneTouchMEA.com to learn if the OneTouch Reveal® Mobile App is available in your country.

Multiple OneTouch Verio Reflect™ Meters can be paired with your compatible wireless device. For example, your compatible wireless device can be paired with a meter at home and another at work. To pair multiple meters, repeat the pairing instructions for each meter. See page 110 for pairing instructions.

Your OneTouch Verio Reflect™ Meter can be paired with multiple compatible wireless devices. To pair multiple compatible wireless devices, repeat the pairing instructions for each compatible wireless device.

Pairing instructions

The BLUETOOTH® feature must be turned on in order to pair the meter with the app. See page 106.

1. Open the OneTouch Reveal® Mobile App and follow instructions to pair the meter with your compatible wireless device

2. Look for "OneTouch" and the last 4 characters of the meter serial number on the compatible wireless device display to correctly identify your meter

The meter screen will display the pairing info to look for on your compatible wireless device.


Take Note

Now use the OneTouch Reveal mobile app and select OneTouch 4444.

OK

Example

3. When prompted by the OneTouch Reveal® Mobile App, the meter will display a six digit PIN number

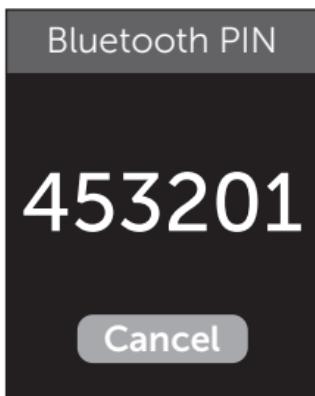

Example PIN

Enter the PIN number into your compatible wireless device using the keypad.

⚠ CAUTION:

Make sure the PIN you enter on your compatible wireless device matches the PIN on your meter display. If a PIN number unexpectedly appears on your meter display, cancel the PIN request by either inserting a test strip to take a test or press the button to cancel.

7 Syncing your meter

4. Wait for your compatible wireless device to indicate that it is paired with your meter

Sending your results to the app

If the BLUETOOTH® feature on the meter is turned on, indicated by the BLUETOOTH® symbol (BT), the meter will automatically send your results to any paired compatible wireless device.

The Sync symbol (C) flashes on the meter display. After syncing, the Sync symbol will disappear and the app will display a list of any new results sent from the meter. You can then share your results with your caregiver or healthcare professional.

⚠ CAUTION:

Always use the current result on your meter for immediate treatment decisions.

NOTE: The compatible wireless device must have the app open and have already been paired to the meter before sending a result. See page 108.

NOTE: If the BLUETOOTH® feature on the meter is **turned off**, or the meter is out of range, the result is not sent to the compatible wireless device. The result is saved in the meter memory with the current date and time, and will be sent to the app the next time you sync. The sent results are also stored in the meter. To sync, the app must be open and running on your compatible wireless device.

To ensure that glucose test results are successfully sent to the app, turn on the BLUETOOTH® feature and check the following:

- The compatible wireless device and meter are both turned on, and the app is running.
- The meter is correctly paired with your compatible wireless device.
- The BLUETOOTH® feature on both devices is running (indicated by ) and the devices are within 8 metres of each other.
- The meter will attempt to transmit results up to 4 hours after a test, even if the meter appears to be off.

7 Syncing your meter

If you are still unable to send results to the compatible wireless device, try un-pairing and re-pairing with the device, moving the meter and compatible wireless device closer together, or relocating them from possible sources of interference. If you continue to have difficulties, please call Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet. Please address any app issues to the app customer service department.

NOTE:

- Inserting a test strip during the transmission will cancel the transfer of all results. The **Apply Blood** screen appears and you can proceed with testing.
- Keep your meter in your care to avoid inaccurate results from being sent to your compatible wireless device.

Downloading results to a computer

Your meter can work with diabetes management software, which provides a visual way to track key factors that affect your blood sugar. To learn more about diabetes management tools available to you, contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.

Connect only to a computer certified to UL 60950-1 (UL).


To transfer meter data, follow the instructions provided with the diabetes management software product to download the results from the meter. You will need a standard micro USB interface cable to connect your OneTouch Verio Reflect™ Meter to a computer to download results (not included).

7 Syncing your meter

Once the command to start the download is sent from the computer to the meter, the **Cable Connected** screen will appear, indicating that the meter is in communication mode.

Do Not insert a test strip while the meter is connected to a computer.

If you are unable to download your results to a computer, please call Customer Service.

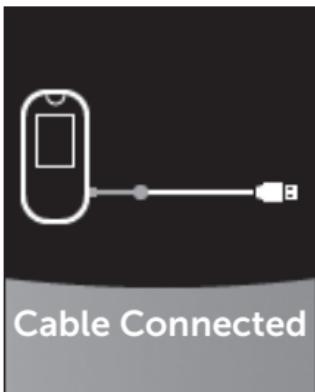

This page left blank intentionally.

8 Adjusting meter settings

Settings Menu overview

After first time setup, you can use the Settings Menu to adjust the meter settings and access additional features.

Time/Date	Set the time and date
Language	Set the language displayed on screens
Range/90d	<ul style="list-style-type: none">• Set Before Meal limits• Set After Meal limits• Set the 90 Day Average Target
Features	<ul style="list-style-type: none">• Turn Basic Mode on or off• Turn Test Tracker on or off• Set Daily Test Goal• Set Range Indicator icon to emoji or text
Messages	<ul style="list-style-type: none">• Turn Pattern Messages on or off• Turn Mentor Tips on or off• Turn Awards on or off
Bluetooth	Turn the BLUETOOTH® feature on or off
Control Test	Initiate a control solution test
Meter Info	View meter serial number, software version and last sub-error code

NOTE: For information on **Bluetooth**, see page 12.
For information on **Control Test**, see page 68.

Navigating to the Settings Menu

From the Main Menu,
press  or  to highlight
Settings and press 

The Settings Menu is displayed.



Adjusting the time and date

You can change the time and date setting whenever needed, such as during travel.

NOTE: Adjusting this setting to a time and date prior to your last blood glucose or control solution test will cause the Results Log to appear out of sequence. The Results Log will display results in the order in which they were taken.

8 Adjusting meter settings

1. From the Settings Menu, press  or  to highlight Time/Date and press 
2. Highlight the time or date and press 
3. Adjust the desired setting and press 

Repeat this step as necessary.

Adjusting the language

1. From the Settings Menu, press  or  to highlight Language and press 

A  appears next to the current language set in the meter.

2. Highlight the desired language, then press 

The Settings Menu will be displayed in the new language you have selected.

Adjusting range limits and 90 Day Average Target


NOTE:

- If you adjust your range limits, your previous Range Indicator icons in the Results Log will not change. Only new results will be affected by your range limit changes.
- Consult with your healthcare professional about the range limits and 90 Day Average Target that are right for you.


8 Adjusting meter settings

1. From the Settings Menu, press or to highlight Range/90d (or Range for Basic Mode) and press

Your current meter settings are displayed.

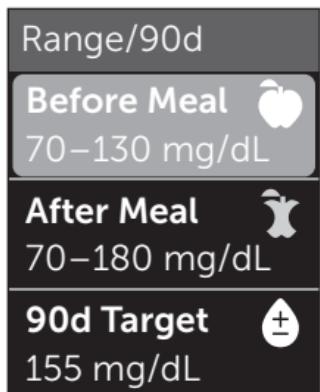


Example - Basic
Mode Off

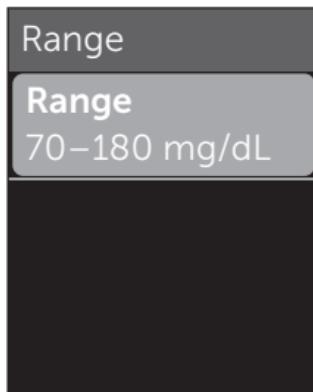

Example - Basic
Mode On

2. Highlight the range setting to be adjusted and press

**3. Adjust the Low Limit
and press **

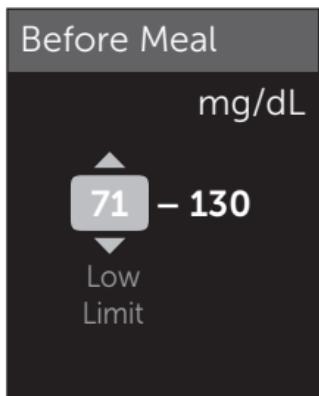

**Example - Before
Meal**

Repeat this step to change the High Limit. See page 27 for information about Before and After Meal low and high limits.

NOTE: In Basic Mode, range limits are general and are not defined by mealtime. General range limits are pre-set.


8 Adjusting meter settings

General range limits

pre-set low limit	70 mg/dL
pre-set high limit	180 mg/dL
low limit must be set within	60 mg/dL and 110 mg/dL
high limit must be set within	90 mg/dL and 300 mg/dL

4. To adjust the 90 Day Average Target, highlight the 90d Target setting and press 

5. Adjust the 90 Day Average Target and press 

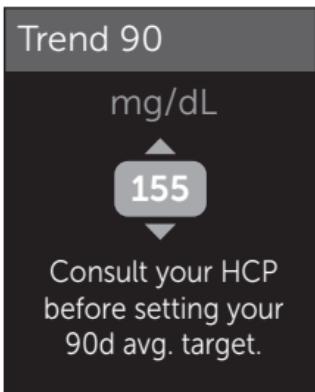

The meter is pre-set with a 90 Day Average Target of 155 mg/dL. The minimum 90 Day Average Target is 99 mg/dL and the maximum is 265 mg/dL.

NOTE: Be sure to talk to your healthcare professional about the 90 Day Average Target that is right for you.

Turning Basic Mode on or off

NOTE: Selecting Basic Mode affects the information the meter displays. See page 37 for more information.

1. From the Settings Menu, press  or  to highlight Features and press 

The Features Menu indicates whether Basic Mode is currently on or off.

2. Highlight Basic Mode and press 

A  appears next to the current mode set in the meter.

8 Adjusting meter settings

3. Highlight On or Off and press

In Basic Mode, Test Tracker and Range Indicator are removed from the Features Menu.

Basic Mode uses pre-set general range limits which are different from the ones set during first time setup. See page 121 for information on adjusting range limits.

Turning the Test Tracker on or off and setting a Daily Test Goal

1. From the Settings Menu, press or to highlight Features and press

The Features Menu indicates the current Test Tracker setting.

NOTE: If Basic Mode is selected the Test Tracker will not appear in the Features Menu.

2. Highlight Test Tracker and press

3. Adjust the Daily Test Goal, then press

The meter is pre-set with a Daily Test Goal of four tests per day. The maximum Daily Test Goal is eight.

NOTE: Be sure to talk to your healthcare professional about the Daily Test Goal that is right for you.

To turn off the Test Tracker, press  or  until **Off** is highlighted, and press .

Set the Range Indicator icon to emoji or text

1. From the Settings Menu, press or to highlight Features and press

Your current meter settings are displayed.

2. Highlight Range Indicator and press

A  appears next to the current setting in the meter.

3. Highlight Emoji or Text, then press

Turning Pattern Messages, Mentor Tips and Awards on or off

NOTE:

- Basic Mode must be turned off for the Messages Menu to appear in Settings. See page 125.
- If Pattern Messages are turned off, you will no longer see Patterns with results, Summary or Results Log screens. The meter will continue to track Pattern information. If you turn Pattern Messages back on again, Pattern icons will be displayed in the Summary and Results Log. See page 83 for more information about Pattern Messages.
- If Mentor Tips are turned off, you will no longer see Mentor Tips with results. See page 90 for more information about Mentor Tips.
- If Awards are turned off, you will no longer see Awards with results or in the Summary view. See page 93 for more information about Awards.

1. From the Settings Menu, press or to highlight Messages and press

Your current meter settings are displayed.

2. Highlight Pattern, Mentor Tips or Awards and press

A appears next to the current setting in the meter.

3. Highlight On or Off and press

View meter info and last sub-error code

The meter serial number, software version, and information about the last meter sub-error are stored in your meter. You can check this information at any time and use it for troubleshooting.

8 Adjusting meter settings

1. From the Settings Menu, press  or  to highlight Meter Info and press 

The meter serial number, software version, and unit of measure are displayed.

Meter Info
Serial Number:
ABC12345
Software:
PUE_01.01.00
Units:
mg/dL

2. Press  to view information about the last meter error

Meter Info
Last Error:
Error X
Error Date:
30 Aug 2018
Sub-Error Code:
XXXX

Pressing  returns you to the Settings Menu.

Example

Example

This page left blank intentionally.

9 Caring for your system

Storing your system

Store your meter, test strips, control solution and other items in your carrying case. Keep in a cool, dry place between 5°C and 30°C. **Do Not** refrigerate. Keep all items away from direct sunlight and heat. Exposure to temperatures outside the storage limits, as well as moisture and humidity, may cause inaccurate readings.

Cleaning and disinfection

Cleaning and disinfection are different and both should be performed. Cleaning is part of your normal care and maintenance and should be performed prior to disinfection, but cleaning does not kill germs. Disinfection is the only way to reduce your exposure to disease. For cleaning information, see page 133, and for disinfection, see page 135.

Cleaning your meter, lancing device and cap

The meter, lancing device and cap should be cleaned whenever they are visibly dirty and before disinfection. Clean your meter at least once per week. For cleaning obtain regular strength liquid dish soap and a soft cloth. Prepare a mild detergent solution by stirring 2.5 mL of regular strength liquid dish soap into 250 mL of water.

- **Do Not** use alcohol or any other solvent.
- **Do Not** allow liquids, dirt, dust, blood or control solution to enter the test strip port or the data port. (See page 20.)
- **Do Not** spray cleaning solution on the meter or immerse it in any liquid.



9 Caring for your system

1. Holding the meter with the test strip port pointed down, use a soft cloth dampened with water and mild detergent to wipe the outside of the meter and lancing device



Be sure to squeeze out any excess liquid before you wipe the meter. Wipe the outside of the cap.

2. Wipe dry with a clean, soft cloth



Disinfecting your meter, lancing device and cap

The meter, lancing device and cap should be disinfected periodically. Clean your meter, lancing device and cap prior to disinfecting. For disinfecting, obtain regular household bleach (*containing a minimum of 5.5% sodium hypochlorite as the active ingredient*)*. Prepare a solution of 1 part household bleach and 9 parts water.

*Follow manufacturer's instruction for handling and storage of bleach.

1. Hold the meter with the test strip port pointed down

Use a soft cloth dampened with this solution to wipe the outside of the meter and lancing device until the surface is damp. Be sure to squeeze out any excess liquid before you wipe the meter.



9 Caring for your system

2. After wiping, cover the surface you are disinfecting with the soft cloth dampened with the bleach solution for 1 minute

Then wipe with a clean, damp, soft cloth.



Wash hands thoroughly with soap and water after handling the meter, lancing device and cap.

If you see signs of wear, please contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.

This page left blank intentionally.

Replacing the batteries

Your OneTouch Verio Reflect™ Meter uses 2x CR2032 lithium batteries. One battery powers the meter only and the other powers the backlight. See page 157 for information on when to change the meter batteries.

If the meter does not turn on, or remains dim after pressing a button, check the batteries.

Do Not change the batteries when connected to a PC.

⚠WARNING: CHEMICAL BURN HAZARD. DO NOT INGEST BATTERY. This product contains a coin/button cell battery. If swallowed, it can quickly cause severe internal burns and can lead to death. Keep new and used batteries away from children. If you think batteries might have been swallowed, seek immediate medical attention.

IMPORTANT: Use only CR2032 lithium batteries with your meter. **Do Not** use rechargeable batteries. Use of an incorrect battery type or a used battery may result in your meter providing fewer tests than normal.

⚠WARNING: Certain batteries may cause leaking, which can damage the meter or cause the batteries to lose power sooner than normal. Replace leaking batteries immediately.

1. Remove the old batteries

Start with the meter turned off. Remove the battery cover by pressing and sliding it downward.



Pull up firmly on the plastic ribbons. The plastic ribbon with the  symbol is for the meter battery, and the plastic ribbon with the  symbol is for the backlight battery. Always change both batteries at the same time.



2. Insert the new batteries

With the "+" side facing up toward you, place each battery in the compartment within the fold of the plastic ribbon.

Push each battery in until it snaps into the battery clasp.



Replace the battery cover by sliding it upwards onto the meter.

If the meter does not power on after you have replaced the batteries, check that the batteries are correctly installed with the "+" side up. If the meter still does not power on, contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.



3. Check your meter settings

You will be prompted to set the time and date whenever you change the batteries. See page 119.

4. Dispose of batteries

Dispose of batteries according to your local environmental regulations.

11 Multi-patient use testing

Healthcare professional information

Only healthcare professionals working in a professional capacity in healthcare settings may perform blood glucose tests on more than one patient using the OneTouch Verio Reflect™ Meter. Always follow recognised procedures for handling objects that are potentially contaminated with human material when using the meter.

CAUTION:

- **Any patient with an infection or suffering from an infectious disease, or is suspected to be, and any patient who is a carrier of a multi-resistant microorganism must be assigned his/her own meter. During this time the meter must not be used to test any other patient.**

- Patients and medical staff are potentially at risk of becoming infected if the same OneTouch Verio Reflect™ Meter is used to test blood glucose in more than one patient. All objects that come into contact with human blood are potential sources of infection.
- When testing multiple patients, it is important to follow your institution's policies and local standards and guidelines that apply to your region.
- Follow your institution's policy and procedure guidelines for when to perform a control solution test.
- Follow your institution's policy for treating hypoglycaemia and hyperglycaemia.
- The patient's hands should be washed with warm water and soap and then dried well prior to testing.
- Results may be falsely low if the patient is severely dehydrated.

11 Multi-patient use testing

- If uncertain about which patient a blood glucose test result belongs to, **retest prior to treatment**.
- **Do Not** use results stored in the Results Log to make immediate treatment decisions. The meter stores up to 750 past results and these may be from multiple patients. Always use the current result for immediate treatment decisions.
- If a single meter is used to test more than one patient, the meter must be **disinfected after each patient**, (see page 135) whether or not blood- or body fluid-contamination is suspected.
- Follow your institution's policy/guidelines for proper disposal of meter and test strips.
- If used for a single patient and a risk of contamination exists, the meter should be disinfected after each use.
- If used for a single patient clean the meter at least once per week.

Personal protective equipment

Healthcare professionals should follow their institution's policy and procedure guidelines for PPE and hand hygiene.

Lancing devices

A single-use, retractable-needle lancing device should be used for each patient skin puncture. Dispose of single-use lancing devices in an appropriate biohazard sharps container immediately after use.

NOTE:

- The single-use lancing device is for getting a capillary whole blood sample only. Not supplied with this kit. If a lancing device and lancets were provided with this kit, dispose of them immediately following your institution's policy and procedure guidelines.

11 Multi-patient use testing

Meter

- It is recommended that the meter is set to Basic Mode for multi-patient use.
- It is important to follow your institution's policies and guidelines before changing any meter settings.
- If possible, a separate OneTouch Verio Reflect™ Meter should be assigned to each patient, and used only to test that patient.
- If the meter is being used with multiple patients, the previous result may belong to a different patient.

Lifespan of meter:

- 9 months in a multi-patient use setting
- During the use of the product, users should properly maintain according to this Owner's Booklet.
- The meter can be disinfected after each test for 9 months and has been tested to 3799 disinfection cycles.

- After 9 months in a multi-patient use setting, as long as you can ensure the basic safety and effectiveness, the meter can continue to be used.
- If the meter stops working or a defect occurs on any part, contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.

Applying blood samples

NOTE: When applying a whole blood sample from the finger, keep the meter pointed down to prevent blood from entering the test strip port.

Follow all instructions for performing a blood glucose test. See page 40.

12 Troubleshooting

Error and other messages

The OneTouch Verio Reflect™ Meter displays messages when there are problems with the test strip, with the meter or when your glucose levels are above 600 mg/dL or below 20 mg/dL. Improper use may cause an inaccurate result without producing an error message.

NOTE: If the meter is on but does not operate (locks up), contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.

What it means

You may have a very low blood glucose level (severe hypoglycaemia), below 20 mg/dL.

What to do

This may require immediate treatment. Although this message could be due to a test error, it is safer to treat first and then do another test. Always treat according to your healthcare professional's recommendations.

⚠ Warning

EXTREME LOW GLUCOSE

Below 20 mg/dL

Treat Low

What it means

You may have a very high blood glucose level (severe hyperglycaemia), over 600 mg/dL.

What to do

Re-test your blood glucose level. If the result is **EXTREME HIGH GLUCOSE** again, obtain and follow instructions from your healthcare professional right away.

⚠ Warning

EXTREME HIGH GLUCOSE

Above 600 mg/dL

12 Troubleshooting

What it means

Meter is too hot (above 44°C) to perform a test.

What to do

Move the meter and test strips to a cooler area. Insert a new test strip when the meter and test strips are within the operating range (6-44°C). If you do not get another **Temperature too high** message, you can proceed with testing.

If this message continues to appear, contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.



Warning

Temperature too high.
Out of operating
range. See Owner's
Booklet.

What it means

Meter is too cold (below 6°C) to perform a test.

What to do

Move the meter and test strips to a warmer area. Insert a new test strip when the meter and test strips are within the operating range (6-44°C). If you do not get another **Temperature too low** message, you can proceed with testing.

If this message continues to appear, contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.



Warning

Temperature too low.
Out of operating
range. See Owner's
Booklet.

Error Screens


If there is a problem with the meter, the following error screens may appear. If you cannot resolve the error with your meter, contact Customer Service. They will refer to the error number, and a sub-error code found in the Meter Info screen, to help troubleshoot the problem. See page 129 for more information on viewing the Meter Info screen.

What it means

There is a problem with the meter.

What to do

Do Not use the meter. Contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.

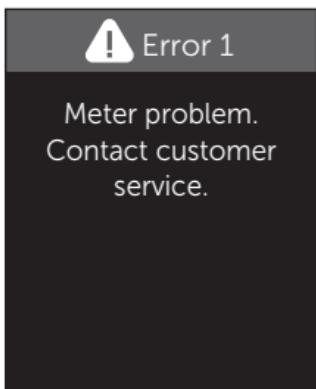

What it means

Error message could be caused by a used test strip, applying blood to the test strip before inserting it into the meter or a problem with the meter or test strip.



Error 2

Meter or strip problem. Retest with a new strip.

What to do

Repeat the test with a new test strip; see page 40 or page 68. If this message continues to appear, contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.

12 Troubleshooting

What it means

The sample was applied before the meter was ready.

What to do

Repeat the test with a new test strip. Apply a blood or control solution sample only after the

Apply Blood or Apply Control

Solution screen appears on the display. If this message continues to appear, contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.



Error 3

Meter was not ready.
Retest with a
new strip.

What it means

One of the following may apply:

- Not enough blood or control solution was applied or more was added after the meter began to count down.
- The test strip may have been damaged or moved during testing.
- The sample was improperly applied.
- There may be a problem with the meter.

What to do

Repeat the test with a new test strip; see page 40 or page 68. If the error message appears again, contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.



Error 4

Strip fill problem.
Retest with a
new strip.

12 Troubleshooting

What it means

The meter has detected a problem with the test strip. Possible cause is test strip damage.

What to do

Repeat the test with a new test strip; see page 40 or page 68. If the error message appears again, contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.



Error 5

Strip problem. Retest with a new strip.

What it means

Battery power is low but there is still enough battery power to perform a test. Test results will still be accurate. The flashing low battery (■) icon will continue to appear until the batteries are replaced.

Be advised



Low Battery.
Replace soon.
2 x CR2032

OK

What to do

Press  to continue, but replace the batteries as soon as possible.

What it means

There is not enough meter battery power to perform a test.

What to do

Replace the batteries immediately.

Warning



Battery Empty.
Replace now.
2 x CR2032

Comparing meter results to laboratory results

Results obtained from the OneTouch Verio Reflect™ Meter and laboratory tests are reported in plasma-equivalent units. However, your meter result may differ from your lab result due to normal variation. A result from your OneTouch Verio Reflect™ Meter is considered accurate when it is within 15 mg/dL of a laboratory method when the glucose concentration is lower than 100 mg/dL and within 15% of a laboratory method when the glucose concentration is 100 mg/dL or higher.

Meter results can be affected by factors that do not affect lab results in the same way. Specific factors that may cause your meter result to vary from your lab result may include:

- You have eaten recently. This can cause a result from fingertip testing to be up to 70 mg/dL higher than a lab test using blood drawn from a vein.¹
- Your haematocrit is above 60% or below 20%.
- You are severely dehydrated.

For additional information, refer to the OneTouch Verio® Test Strip Insert.

¹Sacks, D.B.: "Carbohydrates." Burtis, C.A., and Ashwood E.R. (ed.), *Tietz Textbook of Clinical Chemistry*, Philadelphia: W.B. Saunders Company (1994), 959.

Guidelines for obtaining an accurate meter to lab comparison

Before going to the lab:

- Perform a control solution test to make sure your meter is working properly.
- **Do Not** eat for at least 8 hours before you test your blood.
- Take your meter and testing supplies with you to the lab.

Testing with your OneTouch Verio ReflectTM Meter at the lab:

- Test within 15 minutes of the lab test.
- Use only a fresh, capillary blood sample from your fingertip.
- Follow all instructions in this Owner's Booklet for performing a blood glucose test.

Comparing your meter results to those taken from another meter

Comparing your blood glucose test results taken with this meter to your results taken from a different meter is not recommended. Results may differ between meters and are not a useful measure of whether your meter is working properly.

Technical specifications

Assay method	FAD-GDH (flavin adenine dinucleotide dependent glucose dehydrogenase)
Automatic shutoff	Two minutes after last action
Battery ratings	Two 3.0 V d.c., (2x CR2032 batteries), — 15 mA
Battery type	2 replaceable 3.0 Volt CR 2032 lithium batteries (or equivalent)
BLUETOOTH® technology	Frequency range: 2.4-2.4835 GHz Maximum power: 0.4 mW Operating Range Distance: minimum 8 metres (unobstructed) Operating Channels: 40 Channels Security Encryption: 128-bit AES (Advanced Encryption Standard)
Calibration	Plasma-equivalent

Data port type	Compatible with micro USB 2.0
Memory	750 test results
Operating ranges	6-44°C Relative humidity: Non-condensing 10-90% Altitude: up to 3048 metres Haematocrit: 20-60%
Reported result range	20-600 mg/dL
Sample	Fresh capillary whole blood
Sample volume	0.4 µL
Size	43.0(W) x 101.0(L) x 15.6(T) mm
Test time	Approximately 5 seconds
Unit of measure	mg/dL
Weight	Approximately <53 grams

System accuracy

Diabetes experts have suggested that glucose meters should agree within 15 mg/dL of a laboratory method when the glucose concentration is lower than 100 mg/dL, and within 15% of a laboratory method when the glucose concentration is 100 mg/dL or higher.

Samples from 100 patients were tested using both the OneTouch Verio Reflect™ System and the YSI 2900 Glucose Analyzer laboratory instrument.

System accuracy results for glucose concentrations <100 mg/dL

Percent (and number) of meter results that match the laboratory test

Within ± 5 mg/dL	Within ± 10 mg/dL	Within ± 15 mg/dL
73.7% (137/186)	96.8% (180/186)	100% (186/186)

System accuracy results for glucose concentrations ≥ 100 mg/dL

Percent (and number) of meter results that match the laboratory test

Within $\pm 5\%$	Within $\pm 10\%$	Within $\pm 15\%$
65.2% (270/414)	93.2% (386/414)	99.3% (411/414)

System accuracy results for glucose concentrations between 34.0 mg/dL and 451 mg/dL

Within ± 15 mg/dL or $\pm 15\%$
99.5% (597/600)

NOTE: Where 34.0 mg/dL represents the lowest glucose reference value and 451 mg/dL represents the highest glucose reference value (YSI value).

Regression statistics

Samples were tested in duplicate on each of three test strip lots. Results indicate that the OneTouch Verio Reflect™ System compares well with a laboratory method.

# of Subjects	# of Tests	Slope	Intercept (mg/dL)
100	600	1.02	0.72

95% CI Slope	95% CI Intercept (mg/dL)	Std. Error (S _{y,x}) (mg/dL)	R ²
1.01 to 1.03	-0.62 to 2.07	8.92	0.99

User performance accuracy

A study evaluating glucose values from fingertip capillary blood samples obtained by 313 lay persons showed the following results:

97.1% within ± 15 mg/dL of the medical laboratory values at glucose concentrations below 100 mg/dL, and 96.4% within $\pm 15\%$ of the medical laboratory values at glucose concentrations at or above 100 mg/dL.

96.5% of the total number of samples were within ± 15 mg/dL or $\pm 15\%$ of the medical laboratory values.

Precision

Within Run Precision (300 Venous Blood Samples Tested per Glucose Level)

Data generated using the OneTouch Verio Reflect™ Meter.

Target Glucose (mg/dL)	Mean Glucose (mg/dL)	Standard Deviation (mg/dL)	Coefficient of Variation (%)
40	39.5	1.06	2.68
90	88.7	2.36	2.66
130	128.0	3.82	2.98
200	202.4	5.52	2.73
350	357.6	9.28	2.60

Results show that the greatest variability observed between test strips when tested with blood is 2.36 mg/dL SD or less at glucose levels less than 100 mg/dL, or 2.98% CV or less at glucose levels at 100 mg/dL or above.

Total Precision (600 Control Solution Tests per Glucose Level)

Data generated using the OneTouch Verio Reflect™ Meter.

Glucose Level Ranges (mg/dL)	Mean Glucose (mg/dL)	Standard Deviation (mg/dL)	Coefficient of Variation (%)
Low (25-49)	38.3	1.3	3.3
Mid (96-144)	120.7	2.7	2.2
High (298-403)	348.3	8.2	2.4

Guarantee

LifeScan guarantees that the OneTouch Verio Reflect™ Meter will be free of defects in material and workmanship for three years, valid from the date of purchase. The guarantee extends only to the original purchaser and is not transferable. If the meter stops working or a defect occurs on any part, contact Customer Service. Contact information for Customer Service is listed at the end of this Owner's Booklet.

Electrical and safety standards

This meter complies with CISPR 11:Class B (Radiated Only). Emissions of the energy used are low and not likely to cause interference in nearby electronic equipment. The meter has been tested for immunity to electrostatic discharge as specified in IEC 61326-2-6. This meter complies with immunity to radio frequency interference as specified in IEC 61326-1 and 61326-2-6.

The meter meets the requirements for immunity to electrical interference at the frequency range and test level specified in international standard EN 60601-1-2:2014.

Do Not use the equipment where aerosol sprays are being used, or when oxygen is being administered.

90 Day Average Target 36, 121
After Meal range limits 34, 121
After Meal tag 76
AST 45
Awards 23, 82, 93, 128
Basic Mode 37, 125
Batteries, replacing 138
Battery low icon 7, 157
Before Meal range limits 33, 121
Before Meal tag 76
Blood Sugar Mentor Messages 23, 82
BLUETOOTH® 12, 14, 106, 112
Buttons on meter 21
Cleaning your meter, lancing device and cap 133
Colour Bar 25, 57
ColourSure® Dynamic Range Indicator 25, 33, 57, 59
Comparing meter results to laboratory results 158
Control solution 68
Control solution, discard and expiry dates 67
Control solution, testing 66, 72, 74
Daily Test Goal setting 35, 126
Data port 21
Date setting 31
Dehydration 61, 159
Disinfecting your meter, lancing device and cap 135

Display check	26
Disposal, lancets and test strips.....	66
Downloading results to a computer	115
Emoji	127
Error messages.....	148, 149
Event tagging.....	76
EXTREME HIGH GLUCOSE message	62, 149
EXTREME LOW GLUCOSE message	61, 149
Features	126, 127
Fingertip testing procedure	40
First time setup.....	27
General range limits.....	124
Guarantee.....	170
Healthcare professional	142
Hyperglycaemia	62, 149
Hypoglycaemia	61, 149
Icons.....	4, 8
Infection, reduce the chance	46
Intended use	10
Kit components.....	16
Lancet	44, 47
Lancet, expiry date	46
Lancing device	16, 44, 47, 133
Language.....	29, 120
Meal tagging	76

Mealtime range limits 32, 33, 34
Mentor Tips 23, 82, 90, 128
Messages 83, 128
mg/dL 21, 163
Multi-patient use 142
Pairing 108, 110
Pattern Messages 23, 82, 83, 128
Patterns 83, 84
PIN 111
Plasma calibration 158, 162
Range Indicator icon 4, 7, 21, 38, 57, 59, 78,
102, 118, 121, 126, 127
Range limits 121
Serial number 129
Settings 27, 118, 119, 120, 121, 125, 126, 127, 128
Settings Menu 119
Software, diabetes management 115
Software version 129
Storing your system 132
Sub-error code 129
Summary 94, 100
Symbols 4, 8
Syncing 112
Technical specifications 162
Temperature 40, 41, 74, 132, 150, 151, 163
Test principle 11


Test strip 23, 50, 68
Test strip, applying drop of blood 53, 147
Test strip, discard and expiry dates 42
Test strip, silver prongs 23
Test Tracker 25, 35, 59, 126
Time setting 119
Trend 90 25, 103
Turning meter off 38
Turning meter on 26, 50, 68
Unexpected results 60, 63
Unit of measure 21, 163
Unusual red blood cell count (haematocrit) 159, 163
Welcome screen 26

Notes

Contents covered by one or more of the following U.S. patents: 6,797,150, 6,872,298, 8,163,162, 7,749,371, 8,404,102, 8,449,740, 8,529,751, 8,617,370, 9,575,026, 9,575,027, 10,151,724, 6,602,191, 6,976,958, and 7,156,809. Purchase of this device does not act to grant a use licence under these patents. Such a licence is granted only when the device is used with OneTouch Verio® Test Strip. No test strip supplier other than LifeScan, Inc. is authorised to grant such a licence. The accuracy of results generated with LifeScan meters using test strips manufactured by anyone other than LifeScan has not been evaluated by LifeScan.

LifeScan self-test blood glucose monitoring devices conform to the following EU Directives:

IVDD (98/79/EC):

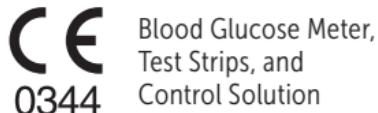

Blood Glucose Meter,
Test Strips, and
Control Solution

MDD (93/42/EEC):



Lancets



Contact LifeScan
Customer Service at:
service@LifeScanMEA.com.



LifeScan Europe GmbH
Gubelstrasse 34
6300 Zug
Switzerland

Meter Made in China



OneTouch®

AW 07187702A 023-923

© 2018–2020 LifeScan

IP Holdings, LLC

Rev. Date: 02/2020

Lifescan™