


NEW

The OneTouch® Delica® Safety Single Use Disposable Lancing Device



Designed for precise and smooth lancing

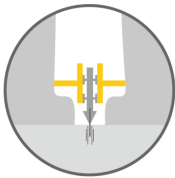

- **Silicone-coated** needles for smooth gliding into the skin
- **Contact-activated** for a consistent sampling experience every time
- **Improved design** for more comfortable testing*
- Cap designed for **easy twist-off** to aid ease of use
- Single use and fully disposable – needle automatically and permanently retracts to **reduce the risk of cross infection**.

Simple 2-step lancing



Twist off and discard lancet cap. **DO NOT PULL.**



Press device firmly against cleaned chosen site to activate.




Dispose of used device in a suitable sharps container.

* Compared to the OneTouch® SureSoft™ single use disposable lancing device.

© 2021 LifeScan IP Holdings, LLC. GL-DLS-2100002

ONETOUCH®

Additional specifications

OneTouch® Delica® Safety 30G Single use disposable lancing device	
Lancing Device	
Needle diameter	30G (0.32 mm)
Penetration depth	1.5 mm
Silicone-coating of needle	Yes
Activation type	Contact activation
Single-use only	Yes
Automatic needle retraction	Yes
Shelf life	5 years from date of manufacturing
Carton variation	<ul style="list-style-type: none">• Box of 200 units• Box of 100 units
Country of manufacture	Japan

Recommend OneTouch® Products today!



service@LifeScanMEA.com



Call OneTouch® Customer Care

U.A.E: 800 4530

Qatar: 800 8512 (Toll Free)

Kuwait: +965 222 76888, ext 220

Bahrain: +973 177 31415

Oman: +968 248 14501

Saudi Arabia: 800 244 0266 (Toll Free)